

# COMPENDIUM OF INTERNATIONAL METHODS OF ANALYSIS FOR SPIRITUOUS BEVERAGES AND ALCOHOLS

OIV-MA-BS-17 Analysis of  $\alpha$ -dicarbonyl compounds by gas chromatography after derivation  
by 1,2-diaminobenzene in spirit drinks of viti-vinicultural origin (Type IV)  
Method OIV-MA-BS-17 : R2010

Type IV method

## Analysis of $\alpha$ -dicarbonyl compounds in spirituous beverages of viti-vinicultural origin by gas chromatography after derivation by 1,2-diaminobenzene

(OIV/OENO 382D/2010)

### 1. Introduction

The principal  $\alpha$ -dicarbonyl compounds found in wine spirits (Figure 1) are: glyoxal, methylglyoxal, diacetyl and pentane-2,3-dione.

Glyoxal	OCH-CHO (ethanedial)
Methylglyoxal	CH <sub>3</sub> -CO-CHO (2-oxopropanal)
Diacetyl	CH <sub>3</sub> -CO-CO-CH (butane-2,3-dione)
Pentane-2,3-dione	CH <sub>3</sub> -CH <sub>2</sub> -CO-CO-CH <sub>3</sub>
Hexane-2,3-dione	CH <sub>3</sub> -CH <sub>2</sub> -CH <sub>2</sub> -CO-CO-CH <sub>3</sub>
<i>Figure 1. The principal <math>\alpha</math>-dicarbonyl compounds of wine (hexane-2,3-dione is not naturally present in wine but it is used as internal standard).</i>	

Dicarbonyl compounds are important because of their sensory impact

### 2. Applicability


This method applies to spirituous beverages of vitivinicultural origin, for a content of carbonyl compounds included between 0.05 mg/L and 20 mg/L.

# COMPENDIUM OF INTERNATIONAL METHODS OF ANALYSIS FOR SPIRITUOUS BEVERAGES AND ALCOHOLS

## OIV-MA-BS-17 Analysis of alpha)-dicarbonyl compounds by gas chromatography after derivation by 1,2-diaminobenzene in spirit drinks of viti-vinicultural origin (Type IV)

### 3. Principle

The method is based on the formation of quinoxaline derivatives from  $\alpha$ -dicarbonyl compounds with 1,2-diaminobenzene (figure 2).

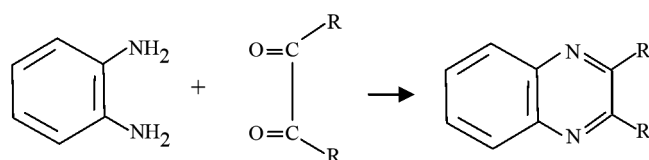

*1,2-Diaminobenzène    Dicarbonyle    Quinoxaline*

*Figure 2 Formation of derivatives.*

The reaction takes place in the spirituous beverage diluted four-fold, pH 8 and after a reaction time of 3 hrs at 60° C. The analysis of the derivatives is then carried out after extraction of the derivatives by dichloromethane and analysis by gas chromatography with detection by mass spectrometry (GC-MS) or using a specific detector of nitrogenous compounds.

### 4. Reagents and products

#### 4.1. Dicarbonyl compounds

4.1.1. Glyoxal (CAS N° 107-22-3) in a 40% solution

4.1.2. Methylglyoxal (CAS N° 78-98-8) in a 40% solution

4.1.3. Diacetyl (CAS N° 431-03-8) > 99 % pure

4.1.4. Pentane-2,3-dione (CAS N° 600-14-6) > 97 % pure

4.1.5. Hexane-2,3-dione (CAS N° 3848-24-6) > 90 % pure

COMPENDIUM OF INTERNATIONAL METHODS OF ANALYSIS FOR SPIRITUOUS BEVERAGES  
AND ALCOHOLS

OIV-MA-BS-17 Analysis of alpha)-diacarbonyl compounds by gas chromatography after derivation  
by 1,2-diaminobenzene in spirit drinks of viti-vinicultural origin (Type IV)

---

**4.2. 1,2-Diaminobenzene (CAS N° 95-54-5) in the forme of powder, > 97 % pure**

**4.3. Water for HPLC (according to standard EN ISO 3696)**

**4.4. Ethanol (CAS N° 64-17-5) pure for HPLC**

**4.5. Sodium Hydroxide (CAS N° 1310-73-2) in a 0,1M solution**

**4.6. Acetic acid (CAS N° 7664-93-9) pure cristallisable**

**4.7. Dichloromethane (CAS N° 75-09-2)**

**4.8. Anhydrous sodium sulphate (n°cas 7757-82-6).**

**4.9. 50% vol. hydroalcoholic solution.**

Mix 50 ml of pure ethanol for HPLC (4.4) with 50 ml of water (4.3)

**4.10. Solution of internal standard hexane-2,3-dione at 2.0 g/l**

Place 40 mg of hexane-2,3-dione (4.2) in a 30 ml flask, dilute in 20 ml of 50% vol. hydroalcoholic solution. (4.9), stir until complete dissolution.

## **5. Apparatus**

**5.1. Gas chromatography with detection by mass spectrometry (GC-MS) or using a special nitrogenised compound detector.**

# COMPENDIUM OF INTERNATIONAL METHODS OF ANALYSIS FOR SPIRITUOUS BEVERAGES AND ALCOHOLS

OIV-MA-BS-17 Analysis of alpha)-diacarbonyl compounds by gas chromatography after derivation

by 1,2-diaminobenzene in spirit drinks of viti-vinicultural origin (Type IV)

**5.1.1. Moderately polar, polyethylene glycol capillary column (such as the Carbowax 20M, BP21) with the following characteristics (as an example) 50 m X 0.32 mm X 0.25 µm.**

**5.1.2. Data acquisition system**

**5.2. pH measuring apparatus**

**5.3. Magnetic stirrer**

**5.4. Mg analytical balance**

**5.5. Oven which can be set to 60°C**

**5.6. Standard laboratory glassware including pipettes, 30-ml screw-cap flasks, and microsyringes**

## **6. Preparation of the sample**

Dilute the spirituous beverage four-fold in water (4.3)

## **7. Procedure**

Place 10 ml of spirituous beverage diluted four-fold (6) in a 30 ml flask

Bring to pH 8 while stirring, with sodium hydroxide 0.1 M (4.5)

Add 5 mg of 1,2-diaminobenzene (4.2)

Add 10 µl of hexane-2,3-dione (internal standard) at 2.0 g/l (4.10)

Close the flask using a screw-cap fitted with a Teflon-faced seal

Stir until the reagent has completely disappeared (5.3)

Place in the oven at 60°C for 3 hrs (5.5)

Cool.

# COMPENDIUM OF INTERNATIONAL METHODS OF ANALYSIS FOR SPIRITUOUS BEVERAGES AND ALCOHOLS

## OIV-MA-BS-17 Analysis of alpha)-diacarbonyl compounds by gas chromatography after derivation by 1,2-diaminobenzene in spirit drinks of viti-vinicultural origin (Type IV)

---

### 7.1. Analysis

#### 7.1.1. Extraction of quinoxalines

- The reactional medium prepared at 7, is brought to pH 2 using H<sub>2</sub>SO<sub>4</sub> 2M (4.6);
- Extract 2 times using 5 ml of dichloromethane (4.7) by magnetic stirring for 5 minutes;
- Elutriate the lower phase each time;
- Mix the two solvent phases;
- Dry on approximately 1 g of anhydrous sodium sulphate (4.8);
- Elutriate.

#### 7.1.2. Chromatographic analysis (given as an example)

- *Detection.* For GC-MS analysis, a Hewlett Packard HP 5890 gas-phase chromatograph was coupled with a chemstation and an HP 5970 mass spectrometer (electron impact 70eV, 2.7 kV),
- *Note:* A specific detector of the nitrogenous compounds can be used.
- *Column.* The column is a BP21 (SGE, 50 m x 0.32 mm x 0.25 µm).
- *Temperatures.* The temperatures of the injector and the detector are 250°C and 280°C, respectively; the oven is kept at 60°C for 1 min., then programmed to increase at a rate of 2°C/min to 220°C, and the final isotherm lasting 20 min.
- *Injection.* The volume injected is 2 µl and the closing time of the injector valves (splitless) is 30s.

#### 7.1.3. Analysis of quinoxalines formed

- *Separation.* The chromatogram of the derivatives by 1,2-diaminobenzene of a wine according to the ion selection method (SIM) is shown in Figure 3.
- *Identification of the peaks.* GC-MS was used to identify the dicarbonyl compounds

# COMPENDIUM OF INTERNATIONAL METHODS OF ANALYSIS FOR SPIRITUOUS BEVERAGES AND ALCOHOLS

## OIV-MA-BS-17 Analysis of alpha)-dicarbonyl compounds by gas chromatography after derivation by 1,2-diaminobenzene in spirit drinks of viti-vinicultural origin (Type IV)

derived from the wine spirit based on the total ionic current method (scan) which is used to obtain the mass spectra of quinoxaline derivatives and to compare them with those memorised in the spectra library; in addition, the retention times were compared with those for pure compounds treated in the same way. Table 1 shows the principal ions of the mass spectra for the obtained dicarbonyl compound derivatives.

- *Proportioning.* The quantitative determination of the dicarbonyl compounds is carried out with the SIM method, by selecting ions  $M/Z = 76, 77, 103, 117, 130, 144, 158$  and  $171$ . The ions  $M/Z = 76$  and  $77$  are used for the quantification and the others as qualifiers, i.e. glyoxal: ions  $M/Z = 103$  and  $130$ , methylglyoxal: ions  $M/Z = 117$  and  $144$ , diacetyl: ions  $M/Z = 117$  and  $158$ , pentan-2,3-dione: ions  $M/Z = 171$  and hexane-2,3-dione: ions  $M/Z = 158$  and  $171$ .

### 7.1.4. Characteristics of the method

- Some internal validation were defined but these are not a formal validation as per the protocol governing the planning, the implementing and the interpreting of the performance studies relating to the analysis methods (OIV 6/2000)
- Repeatability. The repeatability of the GC-MS-SIM method displays variation coefficients ranging between 2 and 5% for the four dicarbonyl compounds;
- Recovery rate. The quantities added to a wine were recovered with with a below 6% deviation from expected results;
- Linearity. Linear correlations were obtained in concentration domains ranging between 0.05 to 20 mg/l.
- Detection limit. The detection limit of most of the derived dicarbonylated products is 0.05 mg/l.

Dicarbonylated compound	Derivative	Mass spectrum (principal ions and abundance)
Glyoxal	Quinoxaline	130 (100), 103 (56.2), 76 (46.8), 50 (20.2)

# COMPENDIUM OF INTERNATIONAL METHODS OF ANALYSIS FOR SPIRITUOUS BEVERAGES AND ALCOHOLS

## OIV-MA-BS-17 Analysis of alpha)-dicarbonyl compounds by gas chromatography after derivation by 1,2-diaminobenzene in spirit drinks of viti-vinicultural origin (Type IV)

Methylglyoxal	2-Methylquinoxaline	144 (100), 117 (77.8), 76 (40.5), 77 (23.3), 50 (21.9), 75 (11.3), 145 (10.3)
Diacetyl	2,3-Diméthylquinoxaline	117 (100), 158 (75.6), 76 (32.3), 77 (23.1), 50 (18.3), 75 (10.4)
Pentane-2,3-dione	2-Ethyl-3-methylquinoxaline	171 (100), 172 (98), 130 (34.1), 75 (33.3), 77(21) 50 (19.4), 144 (19), 143 (14.1), 103 (14)
Hexane-2,3-dione	2,3- Diethylquinoxaline	158 (100), 171 (20.1), 76 (13.7), 77 (12.8), 159 (11.4), 157 (10.8), 50 (8.1)
<i>Table 1. Mass spectra (ion <math>m/Z</math> (intensity of the molecular ion in relation to that of the basic peak) of derivatives of dicarbonyl compounds by 2,3-diaminobenzene</i>		

# COMPENDIUM OF INTERNATIONAL METHODS OF ANALYSIS FOR SPIRITUOUS BEVERAGES AND ALCOHOLS

## OIV-MA-BS-17 Analysis of alpha)-diacarbonyl compounds by gas chromatography after derivation by 1,2-diaminobenzene in spirit drinks of viti-vinicultural origin (Type IV)

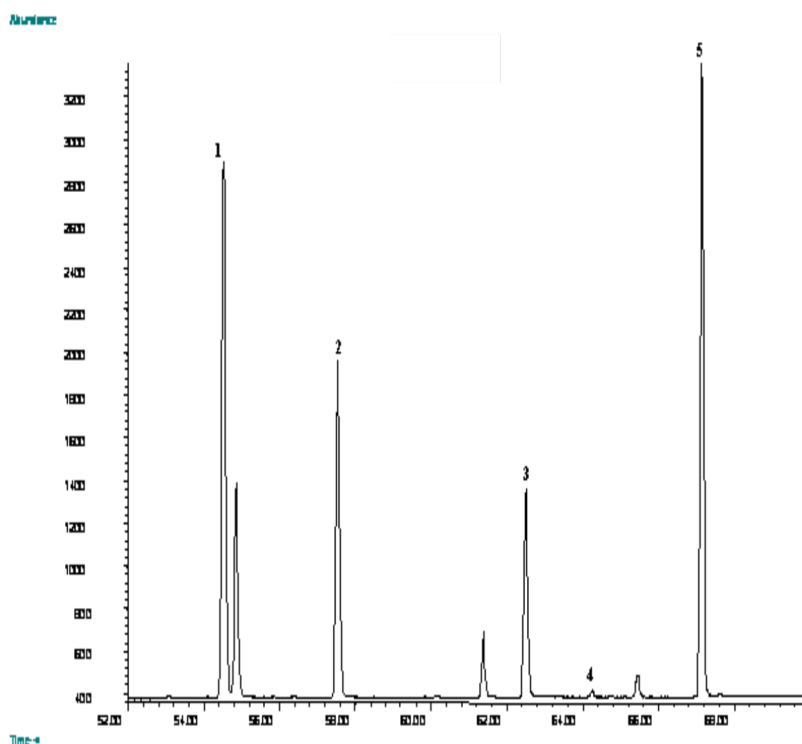


Figure 3. High-performance liquid phase chromatogram of dicarbonyl compounds derived by 1,2-diaminobenzene from wine spirit, detected through mass spectrometry by selecting the ions  $m/z = 76, 77, 103, 117, 130, 131, 144, 158, 160$  and  $171$ . BP21 column,  $50\text{m} \times 0.32\text{mm}$   $0.25\text{ }\mu\text{m}$  oven Temperature  $60\text{ }^{\circ}\text{C}$  for 1 min, then programmed to increase  $2\text{ }^{\circ}\text{C}/\text{min}$  up to  $220\text{ }^{\circ}\text{C}$ . injector Temperature:  $250\text{ }^{\circ}\text{C}$ .

1. Glyoxal
2. Methylglyoxal
3. Diacetyl
4. Pentane-2,3-dione
5. Hexane-2,3-dione (internal standard)
6. Phenylglyoxal (not determined for this method)

## 8. Bibliography

1. Bartowski E.J. and Henschke P.A. The buttery attribute of wine - diacetyl - desirability spoilage and beyond. *Int. j. food microbial.* 96 : 235-252 (2004).



# COMPENDIUM OF INTERNATIONAL METHODS OF ANALYSIS FOR SPIRITUOUS BEVERAGES AND ALCOHOLS

## OIV-MA-BS-17 Analysis of $\alpha$ -dicarbonyl compounds by gas chromatography after derivation by 1,2-diaminobenzene in spirit drinks of viti-vinicultural origin (Type IV)

2. Bednarski W, Jedrychowski L, Hammond E and Nikolov L, A method for determination of  $\alpha$ -dicarbonyl compounds. J Dairy Sci 72:2474-2477 (1989).
3. Leppanen O, Ronkainen P, Koivisto T and Denslow J, A semiautomatic method for the gas chromatographic determination of vicinal diketones in alcoholic beverages. J Inst Brew 85:278- 281 (1979).
4. Martineau B, Acree TE and Henick-Kling T, Effect of wine type on the detection threshold for diacetyl. Food Res Int 28:139-143 (1995).
5. Moree-Testa P and Saint-Jalm Y, Determination of  $\alpha$ -dicarbonyl compounds in cigarette smoke. J Chromatogr 217:197-208 (1981).
6. de Revel G Pripis-Nicolau L. Barbe J.-C; and Bertrand A, The detection of  $\alpha$ -dicarbonyl compounds in wine by formation of quinoxaline derivatives. J Sci. Food Agric.80:102-108 (2000).
7. de Revel G and Bertrand A, Dicarbonyl compounds and their reduction products in wine. Identification of wine aldehydes. Proc 7th Weurman Flavour Research Symp, Zeist, June, pp 353-361 (1994).
8. de Revel G and Bertrand A, A method for the detection of carbonyl compounds in wine: glyoxal and methylglyoxal. J Sci. Food Agric 61:267-272 (1993).
9. Voulgaropoulos A, Soilis T and Andricopoulos N, Fluorimetric determination of diacetyl in wines after condensation with 3,4-diaminoanisole. Am J Enol Vitic 42:73-75 (1991).