

## **RESOLUTION OIV-OENO 406-2011**

### **METHOD FOR THE DETERMINATION OF THE SIZE OF PIECES OF OAK WOOD BY SCREENING**

THE GENERAL ASSEMBLY

In view of Article 2, paragraph 2, no. iv of the Agreement of 3 April 2001 establishing the International Organisation of Vine and Wine,

On the proposal of the Sub-Commission on Methods of Analysis,

DECIDES on the proposal of Commission II "Oenology" to include the following method in the Compendium of international methods of analysis:

#### **Determination of the size of pieces of oak wood by screening**

Type of method: I

#### **1. Introduction**

The use of pieces of oak wood, commonly called chips, to treat wine is authorised provided they comply with the specifications of the Oenological Codex (resolution OENO 3/2005). In particular, the pieces of oak wood used must meet a size requirement, and it is specified that "The dimensions of these particles must be such that at least 95% in weight be retained by the screen of 2 mm mesh (9 mesh)". The following operating procedure provides a method of sampling and then screening that can be used to verify this requirement.

#### **2. Field of application**

The method applies to oak wood test samples of more than 0.5 kg.


#### **3. Principle**

After dividing up the initial test sample, a known quantity of pieces of oak wood (approximately 200g) is placed on a vibrating screen. By weighing the pieces of oak wood remaining on the screen after shaking, it is possible to determine the percentage by weight of particles retained by the screen.

## 4. Equipment

- Standard laboratory equipment.
- Screen of 2 mm mesh (9 mesh), 30 cm in diameter, mounted on a vibrating plate provided with a recovery tray.
- Weighing machine capable of weighing to within 0.1 g.
- Slotted test specimen divider (see figure below as an example).

EN 1482-1:2007 (F)

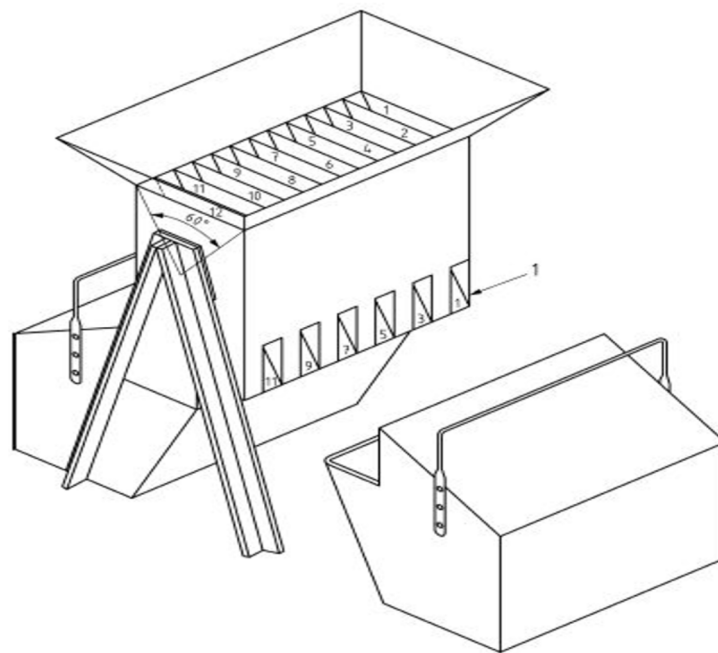

Caption:

1 Alternating sections on either side

Slotted test sample divider (EN 1482-1: 2007)

Scheme proposed as an example

## 5. Division of test sample

When the size of the test sample has to be reduced to obtain “sub-samples” of 200 g which retain a homogeneous nature representative of the initial test sample, a slotted test sample divider can be used which allows random separation of the test sample into 2 parts.

The test sample is poured entirely into the divider in order to separate it into two statistically equivalent parts. Half is put aside, while the other half is again split by means of the chip spreader. This operation is repeated as often as necessary, half being eliminated at each stage with the aim of obtaining 2 “sub-samples” of about 200 g each.

## 6. Operating procedure

- Weigh the empty screen (WES).
- Weigh the empty recovery tray (WET).
- Tare the screen + recovery tray unit and place on it about 200 g of pieces of oak wood weighed to within 0.1 g. Let WOAK be the weight of the pieces of oak wood to be screened.
- Place the unit on the vibrating plate and close the cover with the clamping loops.
- Start up the device and allow it to vibrate for 15 minutes.
- Weigh the screen containing the remaining particles that have not passed through the 2mm meshes (WPS).
- Weigh the recovery tray containing the particles that have passed through the screen (WPT).

A second test is performed in these conditions on the second sub-sample of pieces of oak wood coming from the same initial test sample.

Comment: Weighing of the recovery tray before and after screening (WRT and WPT) serves to verify that there has been no loss of test sample during the operation.

One should have:  $W_{ES} + W_{ET} + W_{OAK} = W_{PS} + W_{PT}$

## 7. Calculation

The percentage (by weight) of particles retained by the screen of 2mm mesh is given by the following formula:

$$\% \text{ of particles retained} = \frac{(W_{PS} - W_{ES}) \times 100}{W_{OAK}}$$

This calculation is performed for each of the 2 sub-samples coming from the initial test sample; the percentage of particles retained corresponds to the mean of the 2 results.

## 8. Bibliography

1. Resolution OENO 3/2005 PIECES OF OAK WOOD
2. EN1482-1 - Fertilizers and liming materials. Sampling and sample preparation. Part 1: Sampling.