

RESOLUTION OENO 16/2003

ACTIVE DRY YEASTS

THE GENERAL ASSEMBLY,

CONSIDERING Article 5, paragraph 4 of the International Convention of the Unification of Methods of Analysis and Appraisal of Wines of 13 October 1954,

UPON THE PROPOSAL of the Sub-commission of Methods of Analysis and Appraisal of Wines,

DECIDES to replace the existing monograph by the following monograph in the International Oenological Codex:

ACTIVE DRY YEASTS (A.D.Y.)

1. OBJECT, ORIGIN AND FIELD OF APPLICATION

Yeasts are used for the inoculation of musts and wine.

The rate of inoculation is at the user's discretion.

- Yeasts used must be isolated from grapes, musts or wine or cultures originating from the combination of these same yeasts (original mother cultures) which must be stored in genetically stable conditions. The obtaining and usage of genetically modified oenological yeast (G.M.Os) require advance authorisation from the competent authority.

2. LABELLING

The following information must be indicated on the label:

- The genus and species name in addition to the reference of the strain(s) attributed by an official body of monitoring micro-organisms or by international institutions, the breeder, the origin and strain breeder and possibly the originator that isolated it.
- Operating instructions or method and the possible reactivation additives recommended by the manufacturer.

- The number of viable cells per gram of powder (CFU as determined in the annex) guaranteed by the manufacturer, loss of viability per month of storage under defined conditions for temperature, humidity and aeration, batch number in addition to the expiration date and storage conditions.
- Where relevant, the indication that the yeasts were obtained through genetic modifications and their modified character.

3. CHARACTERISTICS

Active dry yeast is in the form of round or vermiculated pellets obtained by drying a concentrated yeast culture.

4. TEST TRIAL METHODS AND LIMITS

4.1. Humidity

Measured by the weight loss of 5 g of product dried at 105°C until it reaches a constant weight (about 3 hours).

Maximum level should be less than 8%.

4.2. Heavy metals

Proceed with the determination according to the method in chapter II of the International Oenological Codex.

Content should be less than 10 mg/kg of dry matter, expressed in lead.

4.3. Lead

Proceed with the determination according to the method in chapter II of the International Oenological Codex.

Content should be less than 5 mg/kg of dry matter.

4.4. Mercury

Proceed with the determination according to the method in chapter II of the International Oenological Codex.

Content should be less than 1 mg/kg of dry matter.

4.5. Arsenic

Proceed with the determination according to the method in chapter II of the International Oenological Codex.

Content should be less than 3 mg/kg of dry matter.

4.6. Cadmium

Proceed with the determination according to the method in chapter II of the International Oenological Codex.

Content should be less than 1 mg/kg of dry matter.

4.7. Mycotoxins^[1]

4.8. Viable yeasts

Proceed with counting according to the method in chapter II of the International Oenological Codex (method in annex of this resolution).

Content should be above or equal to 10^{10} CFU/g.

4.9. Yeasts of a different species of the indicated strain

Proceed with counting according to the method in chapter II of the International Oenological Codex (method in annex of this resolution).

Content should be less than 0.01 % of total viable yeasts.

4.10. Moulds

Proceed with counting according to the method in chapter II of the International Oenological Codex (method in annex of this resolution).

Content should be less than 10^3 CFU/g of powder.

4.11. Lactic bacteria

Proceed with counting according to the method in chapter II of the International Oenological Codex (method in annex of this resolution).

Content should be less than 10^4 CFU/g.

4.12. Acetic bacteria

Proceed with counting according to the method in chapter II of the International Oenological Codex (method in annex of this resolution).

Content should be less than 10^3 CFU/g.

4.13. Salmonella

Proceed with counting according to the method in chapter II of the International Oenological Codex (method in annex of this resolution).

Absence should be checked on a 25 g sample.

4.14. *Pseudomonas aeruginosa*^[2]

4.15. *Escherichia coli*

Proceed with counting according to the method in chapter II of the International Oenological Codex (method in annex of this resolution).

Absence should be checked on a 1 g sample.

4.16. Staphylococci

Proceed with counting according to the method in chapter II of the International Oenological Codex (method in annex of this resolution).

Absence should be checked on a 1 g sample.

4.17. Coliforms

Proceed with counting according to the method in chapter II of the International Oenological Codex (method in annex of this resolution).

Number should be less than 10 CFU/g

5. ADDITIVES

They must be in conformity with regulations in force.

6. STORAGE CONDITIONS

Storage should not be in open packaging and/or at temperatures above 10 C.

Storage conditions differ according to preparation and packaging methods.

Always refer to manufacturer's recommendations.

METHODS OF MICROBIOLOGICAL ANALYSIS

(to appear in chapter II of the International Oenological Codex)

1. Preliminary rehydration of active dry yeasts (ADY)

- weigh 1 g of ADY under sterile conditions;
- add 100ml of sterile water at room temperature (20°C) under sterile conditions;
- slowly homogenise using a rod and a magnetic stirrer for 5 min;
- stop stirring and allow to stand for 20 minutes at a temperature of 25°C-30°C;
- homogenise again at room temperature for 5 minutes;
- take 10 ml under sterile conditions and then proceed with micro-biological controls on the homogenised reference solution.

2. Determine number of viable yeasts

2.1. YM agar medium (MALT WICKERHAM)

Composition:

| | |
|----------------------|---------|
| Bacteriological agar | 15 g |
| Yeast extract | 3 g |
| Malt extract | 3 g |
| Peptone | 5 g |
| Glucose | 10 g |
| Water q.s.f. | 1000 ml |

Prior to use, the medium is autoclaved at 120°C for 20 minutes.

After inoculation, the dishes are incubated at 25°C under anaerobic conditions for 48 to 72 hours.

Count the number of CFU and refer to the weight of the dry matter.

2.2. YMS agar medium

Composition:

| | |
|------------------|---------|
| Agar | 20 g |
| Glucose | 20 g |
| Yeast extract | 5 g |
| Malt extract | 3 g |
| Peptone | 2 g |
| Malic acid | 4 g |
| Grape juice | 100 ml |
| Vitamin complex* | 1 % |
| Water q.s.f. | 1000 ml |

Prior to use, the medium is autoclaved at 120°C for 20 minutes.

After inoculation, the dishes are incubated at 25°C under anaerobic conditions for 48 to 72 hours.

Count the number of CFU and refer to the weight of the dry matter.

* Vitamin complex (inositol 25 mg, biotin 0.02 mg, Ca pantothenate 4 mg, folic acid 0.002 mg, nicotinamide 4 mg, paraminobenzoic acid 2 mg, pyridoxine hydrochloride 4 mg, riboflavin 2 mg, thiamine 10 mg, water q.s.f. 1000)

2.3. OGA medium

Composition:

| | |
|----------------------------------|---------|
| Autolytic yeast extract | 5 g |
| Glucose | 20 g |
| Bacteriological agar agar | 15 g |
| Water q.s.f. | 1000 ml |
| Autoclavage at 120°C for 20 min. | |

After inoculation, aerobiosis incubation at 25°C for 48 to 72 hours.
Count the number of CFU and refer to the weight of the dry matter.

3. Counting of yeasts of a different species of the *Saccharomyces* strain according to the lysine test

Lysine test

The yeasts are cultivated in the medium with lysine whose composition is the following:

| | |
|----------------------------|---------|
| Agar | 20 g |
| L-lysine monohydrochloride | 5 g |
| Glucose | 1 g |
| Bromocresol purple | 0.015 g |
| Water q.s.f. | 1000 ml |

Adjust pH 6.8 \pm 0.2

Autoclave at 120°C for 20 min.

After inoculation, the dishes are incubated at 25°C for 48 to 72 hours.
Count the number of CFU and refer to the weight of the dry matter.

4. Counting mould

Czapeck-Dox/s gelose medium

Composition:

| | |
|---------------------------------|--------------|
| Agar agar | <u>15 g</u> |
| Saccharose | <u>30 g</u> |
| NaNO ₃ | <u>3 g</u> |
| K ₂ HPO ₄ | <u>1 g</u> |
| MgSO ₄ | <u>0.5 g</u> |
| KCl | <u>0.5 g</u> |
| FeSO ₄ | 0.01 g |
| Potassium sorbate | 0.4 g |
| Water q.s.f. | 1000 ml |

Adjust pH 3.5

Sterilisation at 120°C for 20 min.

Add directly to the Petri dish 0.1 ml a 0.25% penicillin solution in pure alcohol.

Aerobic incubation at 20°C for 10 days.

5. Counting lactic bacteria

5.1. MTB/s agar medium

Composition:

| | |
|---------|------|
| Glucose | 15 g |
|---------|------|

| | |
|---|---------|
| Peptone | 8 g |
| Yeast extracts | 5 g |
| Casein hydrolysate | 1 g |
| Tomato juice | 20 ml |
| Na acetate | 3 g |
| NH ₄ citrate | 2 g |
| Malic acid | 6 g |
| Mg sulphate | 0.2 g |
| Mn sulphate | 0.035 g |
| Tween 80 1 ml | |
| TC minimal Eagl vitamin | 10 ml |
| after sterilisation adjust to pH 5.0 and add | |
| Agar | 2 % |
| Water q.s.f. | 1000 ml |

Potassium sorbate (400 mg/l in liquid medium) or

Add directly to the Petri dish 0.2 ml of pimarinic hydroalcoholic solution at 25% m/v

Sterilisation at 120°C for 20 minutes

Anaerobic incubation to contrast moulds at 25°C for 8 to 10 days.

5.2. Milieu Man, Rogosa and Sharpe (MRS)

The bacteria are cultivated in a MRS medium (Man, Rogosa, Sharpe 1960) and the composition is as follows:

| | |
|---------------------------------|---------|
| Agar agar | 15 g |
| Bacto-peptone | 10 g |
| Meat extract | 10 g |
| Yeast extract | 5 g |
| Sodium acetate | 5 g |
| K ₂ HPO ₄ | 2 g |
| Trisodium citrate | 2 g |
| MgSO ₄ at 100 mg | 2.5 ml |
| MnSO ₄ at 20 mg | 2 ml |
| Tween 80 1 ml | |
| DL malic acid | 5 g |
| Concentrated tomato juice* | 20 ml |
| Glucose | 20 g |
| Adjust (HCl or NaOH) | pH 4.8 |
| Distilled water q.s.f. | 1000 ml |

Autoclave at 120°C for 20 min.

Potassium sorbate (400 mg/l in liquid medium) or

Add directly to the Petri dish 0.2 ml of pimarinic hydroalcoholic solution at 25% m/v.

Anaerobic incubation at 25°C for 8 to 10 days.

*tomato juice is used to improve lactic bacterial growth.

preparation: take canned tomato juice containing at least 7 g/l of NaCl (maxi 9 g/l)

centrifuge at 4000 g for 20 min;

gather the clear juice and filter through paper filter;
autoclave at 110°C for 20 min.

6. Count of acetic bacteria

6.1. Act/s agar

Composition:

| | |
|---------------------------|--------|
| Bacteriological agar agar | 20 g |
| Yeast extract | 5 g |
| Casein amino acids | 5 g |
| Glucose | 10 g |
| adjust to | pH 4.5 |

Water q.s.f. 1000 ml

Sterilisation

Aerobic incubation at 25°C for 7 days.

Potassium sorbate (400 mg/l in liquid medium) or

Add directly to Perti dish 0.2 ml of pimarcine hydroalcoholic solution at 25% m/v.

7. Count of Salmonella

7.1. Principle

The sample undergoes a pre-enrichment phase in peptoned buffered water for 16 to 20 hours at 37°C. Then the aliquot part of this mixture is inoculated for culture. This contains a specific medium and 2 special tubes (made up of 2 parts) and is incubated 24 hours at 41°C. *Salmonella* migrates from the bottom (selective medium) to the top part of the tube (indicator medium). The presence of *Salmonella* is indicated by a change in colour of this solution.

7.2. Apparatus and analytical conditions

Preparation for culture is carried out in the sterile zone ensured by the Bunsen burner. The soiled material is submitted for destruction by autoclave for 1 hour at 120°C or by total immersion in a bleaching agent for at least 18 hours (See cleaning procedure).

Sterile glass test tube in 125 ml

sterile stomacher bag

Closing barrette

Stomacher

Sterile glass tubes 16x160 mm.

Cottoned glass test tubes 20x220

2 ml sterile plastic pipettes graduated by 0.1 ml

10 ml sterile plastic pipettes graduated by 0.1 ml

Tube shaker

Method for culture to be rehydrated.

2 ml sterile needle with plastic sterile syringe.

Tweezer forceps

Wrench for unscrewing tubes A and B for culture method

Clean glass slide

Sterile cottoned Pasteur pipettes

Wire hoop

Oven at 41°C ± 1°C

Oven at 37°C ± 1°C

Bunsen burner

7.3. Reagents

Sterile peptoned water (SPW)

Sterile distilled water (SDW)

Sterile 500 ml sealed flask filled with 125 ml of SPW

Sterile 500 ml sealed flask filled with 225 ml of SPW

Special medium for *Salmonella*: SRTEM

Novobiocin disk (1.8 mg of novobiocin)

Hektoën agar agar (see DOMIC-08)

API 20E gallery

Agar agar tubes TSAYE inclined

Sterile NaCl at 8.5 g/l solution

Anti-*Salmonella* serum

7.4. Procedure

7.4.1. Preparation of reference suspension

This differs according to nature of products and dilution rate.

Add a test portion of 25 grams or millilitres of the product in a stomacher bag to a nine fold greater amount of peptoned water.

Close the bag by heat welding or using a barrette.

Grind in a stomacher for 1 minute.

7.4.1.1. Pre-enrichment phase in a non selective liquid medium:

Incubate the reference suspension for 16 to 20 hours at $37^{\circ}\text{C} \pm 1^{\circ}\text{C}$.

7.4.1.2. Enrichment in selected liquid mediums

Preparation of culture measures

- unscrew the lid of the culture container;
- add SDW up to line 1 as marked on the container.

Note: The base of tubes A and B must be located under water level.

- adjust the needle to the syringe and check that the syringe plunger is pushed in (absence of air);
- vertically introduce the needle to the syringe in the rubber disc in the centre of the stopper in tube A (blue stopper). Check that the needle is visible under the stopper;
- carefully withdraw the syringe up until the liquid reaches line 3 on the container.

Note: Do not draw up liquid into the syringe.

This operation should take approximately 5 seconds.

- Repeat this operation with tube B (red stopper);
- Close the stopper from the culture container tightly;
- Press the side of the recipient on a tube shaker and maintain at least 5 seconds.

Note: the liquid in tubes A and B must be shaken vigorously.

- Let the culture at least 5 minutes;
- Unscrew the culture container's stopper and pour in the SRTEM medium until the level reaches line 2 as marked on the container;
- Add a novobiocin disc using a tweezer forceps;
- Remove the stoppers from tubes A (blue) and B (red) using a wrench, then dispose of the stoppers.

Note: avoid touching the tubes and the inside wall of the container.

Inoculation of culture container

- Homogenise the pre-enriched culture;
- Identify the culture container. Write down the analysis number on the lid.
- Unscrew the lid.
- Using a 2 ml pipette introduce 1 ml of pre-enriched culture in the culture container.
- Tighten the lid on the culture container.
- Write down the incubation time and date.
- Incubate 24 hours \pm 30 min at 41°C \pm 1°C in a strictly vertical position.

7.4.2. Reading and interpretation

This is carried out by observing the top part of tubes A and B through the container walls.

The possible presence of *Salmonella* is characterised by modifications in indicator medium colour located in one or both of the top parts of the tubes:

| rEaction | TUBE A | tube B |
|-------------------|--------------------------------|---------------------------------------|
| Positive : | All degrees of black colouring | All degrees of red or black colouring |
| Negative : | Absence of black colouring | Absence of red or black colouring |

Tubes showing a positive reaction are subjected to selective agar isolation.

- Dry boxes of Hektoën agar in an incubator at 46°C \pm 1°C until the drops on the surface of the medium disappear completely (lid removed and agar surface facing down).
- Take a wire hoop from the positive middle indicator and inoculate it into 5 ml of SPW, in a 16x160 mm sterile glass tube in order to dilute the culture.
- proceed as such with each positive tube.
- identify the dish and write down on the lid the number of the analysis and the letter of the tube being confirmed.
- homogenise the culture and take a wire hoop.
- isolate the Hektoën agar on the surface to enable the development of isolated colonies.
- incubate 24 hours at 37°C \pm 1°C.
- Select at least 2 isolated colonies considered to be typical.

7.4.3. Confirmation

7.4.3.1. Biochemical tests

- identify the different colonies by using specific miniaturised galleries (API 20^E gallery) by referring to the recommendations of the manufacturer.
- incubate 24 hours at 37°C \pm 1°C.
- At the same time inoculate: an agar to confirm the purity of the strain.
- 1 agar TSAYE inclined for serological typing.

- incubate 24 hours at 37°C \pm 1°C.
- Read the API20E gallery following the manufacturer's indications.
- compare the profile obtained to the standard profiles given by the manufacturer.
- Store TSAYE agar in the refrigerator until utilisation.

7.4.3.2. Serological tests:

Tests are conducted if the strain profile corresponds to *Salmonella* following the recommendations defined by the manufacturer from cultures obtained on agar and after eliminating self-agglutinating strains.

Elimination of self-agglutinating strains:

- Place a drop of 8.5 g/l saline solution on a perfectly clean glass slide.
- Disperse a little bit of the culture taken from the nutritive agar to obtain a homogeneous and cloudy solution using a Pasteur pipette.
- Oscillate the slide for 30 to 60 seconds.
- On a black background using a magnifying glass: if any observation reveals more or less distinct clusters, the strain is considered as being self-agglutinating and should not be subjected to serological typing.


7.5. Results

According to the results based on the interpretation of biochemical and serological testing, the results are expressed as follows:

- Presence of *Salmonella* in m number of grams or ml of product.
- Absence of *Salmonella* in m number of grams or ml of product.

Diagram of procedures

Diagram of procedures

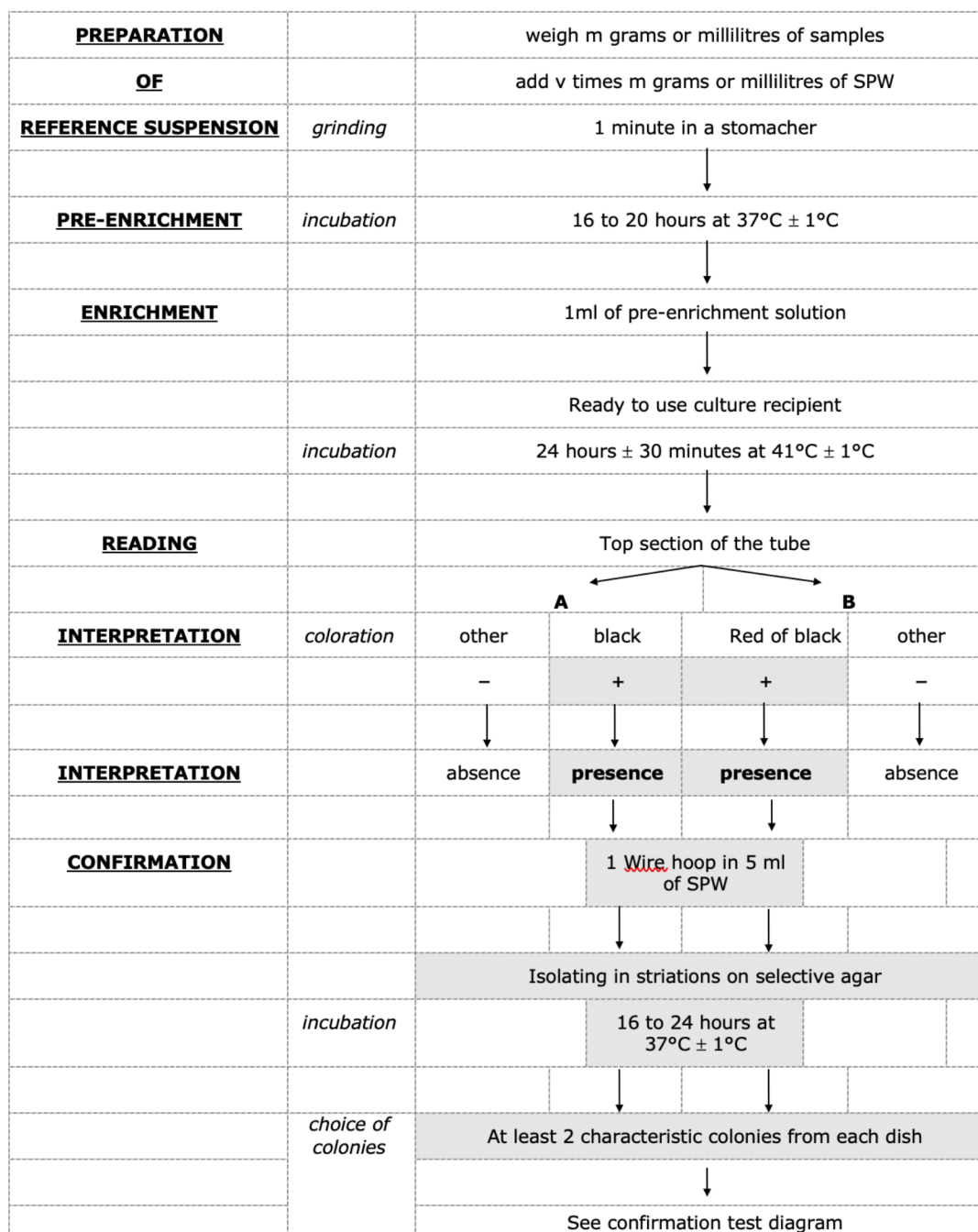

Test confirmation diagram

Test confirmation diagram

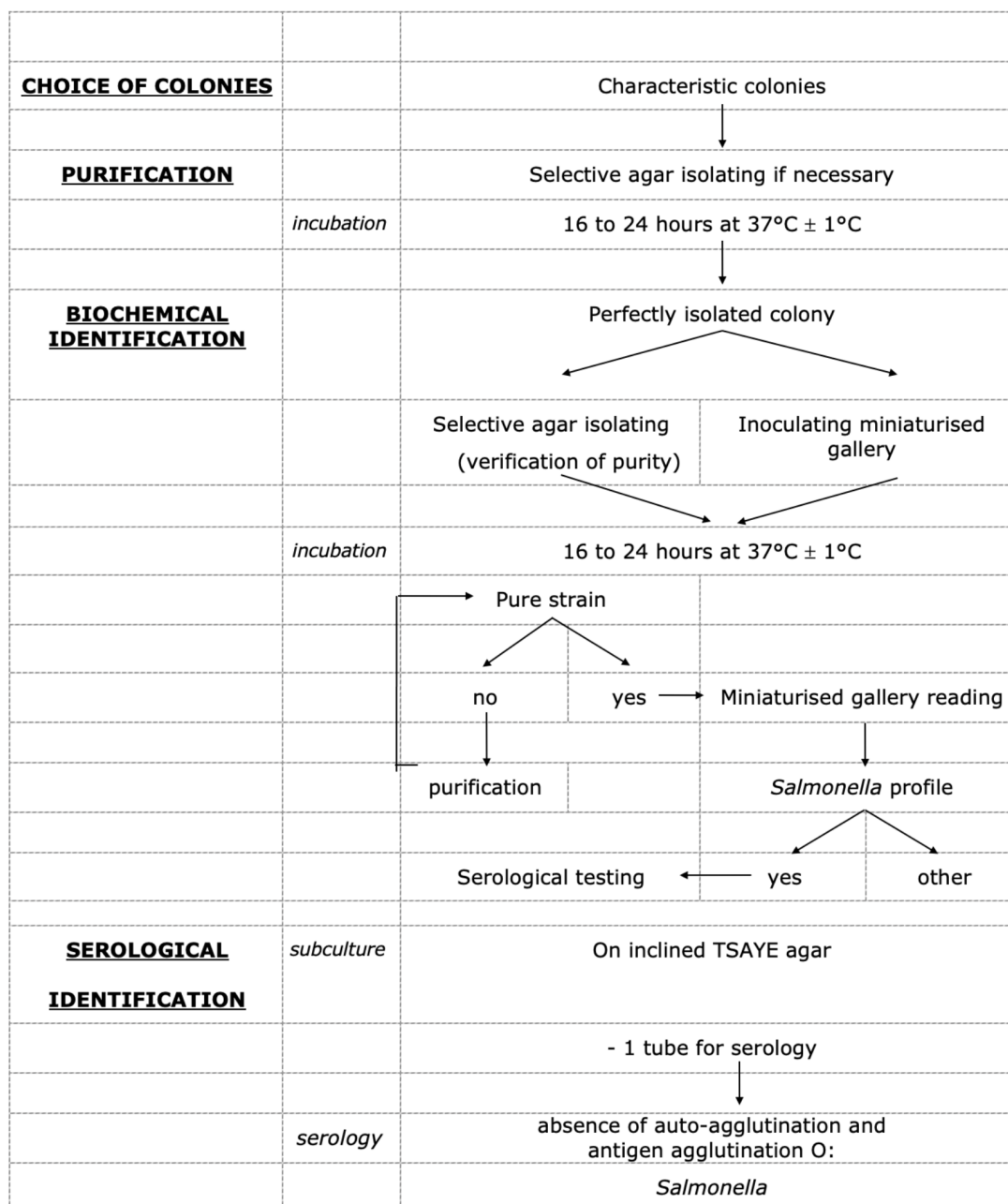


Diagram of biochemical and serological interpretations

| Biochemical reactions | Self-agglutination | Serological reactions | Interpretation |
|-----------------------|--------------------|-----------------------|---|
| typical | no | “O” positive antigen | <i>Salmonella</i> |
| typical | no | Negative reactions | Sent to an authorised centre |
| typical | yes | Not carried out | for determination of the serological type |

8. Count of *Escherichia coli* by the counting of colonies obtained at 44°C

8.1. Principle

Inoculating rapid E. *coli* agar in depth is carried out in a Petri dish for each of the dilutions chosen. Following a 24 hour incubation at 44°C, all characteristic colonies which appear are counted.

8.2. Apparatus and analytical conditions

Cultures are carried out in a sterility zone ensured by the usage of a Bunsen burner.

Plastic sterile Petri dishes with a diameter of 90 millimetres

Sterile 16x60 cottoned glass test tubes

Tube holder

2 ml plastic sterile pipettes with 0.1 ml graduations

Water bath at 100°C ± 2°C

Water bath at 47°C ± 2°C

Tube shaker

Oven at 44°C ± 1°C

Bunsen burner

Colony counter

8.3. Reagents

Sterile diluent for decimal dilutions: tryptone salt (TS)

16x160 pre-filled sterile tubes with 9ml of sterile TS

Rapid *E.coli* cooling agar (R.E.C)

8.4. Procedure

8.4.1. Agar agar medium

- Melt R.EC agar in a boiling water bath. Avoid overheating.
- Never use a culture medium above 50°C.
- For immediate usage, keep agar in the water bath at 47°C ± 2°C.
- Do not cool over 8 hours.
- For a deferred usage maintain the cooling agar in an oven at 55°C ± 1°C.
- The melted culture medium not used within 8 hours will not re-solidify for another usage.

8.4.2. Culture

- Homogenise each dilution before inoculation in Petri dishes and before carrying out decimal dilutions.
- Transfer 1 ml from the reference solution and/or the retained decimal dilutions in the respective Petri dishes. Change the pipette after each dilution.
- Introduce at least 20 minutes after inoculum, 15 to 20 ml of R.EC maintained in the water bath at 47°C ± 2°C.
- Slowly homogenise by shaking.
- Let solidify on the bench (lid up).
- Pour 4 to 5 ml of R.EC maintained at 47°C ± 2°C.
- Let solidify on a bench (lid up).
- Return the dishes and incubate in an oven 24 hours ± 2 hours at 44°C ± 1°C.

8.4.3. Count

Dishes containing between 15 and 150 characteristic colonies of two successive solutions are retained for counting.

If the dish inoculated with 1 ml of first dilution contains characteristic colonies and fewer than 15, it will be retained for counting.

Characteristic colonies are counted using a counter or are counted manually after 24 hours ± 2 hours of incubation.

8.5. Results

8.5.1. General case

The dishes contain between 15 and 150 characteristic colonies for two successive dilutions.

8.5.1.1. Method of calculation

The two dishes retained have between 15 and 150 characteristic colonies. The number N of counted micro-organisms at 44.5°C per millilitre or by gram of product is obtained by calculating the weighted mean on 2 dishes retained.

$$N = \frac{\sum c}{1,1d}$$

$\sum c$: sum of characteristics counted on 2 dishes retained

d : rate of dilution corresponding to first dilution

8.5.1.2. Expression of results

- Round off the number N to 2 significant digits
- Express to the tenth power

ex.: $1.6 \cdot 10^3$ / g or ml

8.5.2. Estimation of small numbers

If the dish inoculated with 1 ml of the 1st retained solution for analysis contains at least 15 characteristic colonies, express the result as follows:

$$N = c \frac{1}{d}$$

c : sum of characteristic colonies counted

d : rate of dilution

If the dish inoculated with 1 ml of the 1st retained solution for analysis does not contain any colonies, express the result as follows:

$$N = < 1 \frac{1}{d} \text{ micro-organism per g or ml}$$

d : rate of dilution

9. Count of Staphylococci with a positive coagulase by the counting and confirmation of colonies obtained at 37°C

9.1. Principle

Decimal dilutions and inoculation on the surface of 1 Baird Parker agar drawn previously in a Petri dish with each of the dilutions retained, are carried out simultaneously from the sample (liquid product) or from the reference solution (other products).

After an incubation of 48 hours at 37°C the characteristic and/or non characteristic colonies are counted and then confirmed by the coagulase test.

9.2. Apparatus and analytical conditions

A. Cultures are carried out in a sterility zone ensured by the usage of a Bunsen burner.

- Sterile glass 16x160 cottoned test tubes
- Sterile plastic precipitating tubes with plastic stoppers
- Tube holder
- 2 ml plastic sterile pipettes with 0.1 ml graduations
- Sterile plastic spreader

- Sterile Pasteur pipettes
- Tube shaker
- Incubate at $37^{\circ}\text{C} \pm 1^{\circ}\text{C}$
- Bunsen burner
- Colony counter

9.2.1. Reagents

- Sterile diluent for tryptone salt (TS) decimal dilutions.
- 16x160 sterile tubes pre-filled with 9ml of sterile TS.
- Baird Parker agar pre-poured in a Petri dish.
- Tubes pre-filled with 5ml brain heart bouillon (sterile).
- Plasma of lyophilised rabbit rehydrated at the time of use.

9.2.2. Procedure

9.2.2.1. Culture

- Dry the agar plates in an incubator at $46^{\circ}\text{C} \pm 1^{\circ}\text{C}$ until the droplets on the surface of the environment have completely disappeared (cover is removed and the agar surface is turned downwards).
- homogenise each dilution prior to inoculation of the surface of agar plate surface before carrying out decimal dilutions.
- Place 0.1 ml of reference solution and/or the retained decimal dilutions on the agar surface while changing the pipette after each dilution.
- Carefully spread the inoculum as quickly as possible using a spreader without touching the edges of the plate.
- Leave the plates with the lids closed for 15 minutes at room temperature.
- Incubate 48 hours \pm 2 hours at $37^{\circ}\text{C} \pm 1^{\circ}\text{C}$

9.2.2.2. Counting

Dishes containing less than 150 characteristic and/or non-characteristic colonies on two successive dilutions are retained, but one of them must include at least 15 colonies. The characteristic and/or non-characteristic colonies are counted either manually or by using a counter.

Characteristic colonies after 48 hours \pm 2 hours of incubation:

- Black or grey, shiny or convex with at least a 1 mm in diameter and a maximum of 2.5 mm in diameter outlined with lightening and precipitation halos.

Non-characteristic colonies after 48 hours \pm 2 hours of incubation:

- Black and shiny with or without a white edge with lightening or precipitation halos absent or barely visible.
- Grey without light zones.

9.2.2.3. Confirmation

Take 3 characteristic colonies or 3 colonies of each type (characteristic or non-characteristic) and submit them to the coagulase test.

Coagulase test:

a. Bouillon culture:

- Take part of the selected colony using a Pasteur pipette sterilised with the Bunsen burner flame and inoculate into a brain heart bouillon.
- Repeat this manipulation for other selected colonies.
- Identify the tubes by sample number and its dilution with a blue marker for characteristic colonies and a green marker for non-characteristic colonies.
- Incubate at $37^{\circ}\text{C} \pm 1^{\circ}\text{C}$ for 20 to 24 hours \square 2H.

b. Testing for free coagulase:

- Add 0.5 ml of culture obtained in brain heart bouillon to 0.5 ml of rehydrated rabbit plasma in a sterile precipitating tube and identify as follows.

- Repeat this procedure for each bouillon culture.
- Incubate 4 to 6 hours at $37^{\circ}\text{C} \pm 1^{\circ}\text{C}$.
- Check for the presence of coagulum or examine the tube after 24 hours \square 2 hours of incubation.

9.2.3. Results

Coagulase is considered positive when it occupies $\frac{3}{4}$ of the initial volume of the liquid.

9.2.3.1. General case

The plates contain a maximum of 150 characteristic and/or non-characteristic colonies.

Calculation procedure:

- Number of Staphylococci with positive coagulase for each plate: a

$$a = \frac{b^c}{A^c} \times c^c + \frac{b^{nc}}{A^{nc}} \times c^{nc}$$

| | |
|----------|---|
| A^c | is the number of spotted characteristic colonies |
| A^{nc} | is the number of spotted non-characteristic colonies; |
| b^c | is the number of characteristic colonies of positive Staphylococci coagulase; |
| b^{nc} | is the number of non-characteristic colonies of positive coagulase Staphylococci |
| c^c | Is the total number of characteristic colonies of positive coagulase Staphylococci for the plate retained; |
| c^{nc} | Is the total number of non-characteristic colonies of positive coagulase Staphylococci positive for the plate retained. |

Round off the number to the nearest whole number.

- Number of positive coagulase Staphylococci in trials: N

The weighted average, calculated as follows from two successive retained solutions:

$$N = \frac{\sum a}{1,1 \times F} \times 10 \text{ positive coagulase Staphylococci by g or ml}$$

$\sum a$: sum of positive coagulase Staphylococci colonies identified on 2 retained plates

F : rate of dilution corresponding to the 1st retained dilution.

Expression of results:

- round off the number N to the two largest whole digits
- express to the tenth power

| ex.: | Amount obtained | Amount rounded off | Result |
|------|-----------------|--------------------|------------------|
| | 36364 | 36000 | $3.6 \cdot 10^4$ |

9.2.3.2. Estimation of small numbers:

If the plate inoculated with 0.1 ml of the first dilution retained for analysis contains less than 15 colonies, the result will be expressed as follows:

$$N = a \frac{1}{d} \times 10 \text{ positive coagulase Staphylococci per g or ml}$$

a : number of positive coagulase Staphylococci identified.

d : rate of dilution for the first dilution retained for analyse.

If the dish inoculated with 0.1 ml of the first dilution retained for analysis contains no positive coagulase Staphylococci the result shall be expressed as follows:

$$N < \frac{1}{d} \times 10 \text{ no positive coagulase Staphylococci per g or ml}$$

d : Rate of dilution from the first retained dilution for analysis.

10. Coliform count by counting colonies obtained at 30°C

10.1. Principle

A. Inoculation in deeply in crystal violet to neutral red (VRBL) lactose bile agar was carried out in Petri dishes for each of the dilutions retained. After incubation for 24 hours at 30°C, the characteristic colonies were counted.

10.2. Apparatus and analytical conditions

Cultures are carried out in a sterile environment as ensured by a Bunsen burner.

- Plastic sterile Petri dishes with a diameter of 90 millimetres
- Sterile glass 16 x 160 cottoned tubes
- Tube holder
- 2 ml plastic sterile pipettes graduated at 0.1 ml
- Water bath at 47°C \pm 2°C.
- Tube shaker
- Incubate at 30°C \pm 1°C
- Incubate at 55°C \pm 1°C
- Bunsen burner
- Colony counter

10.3. Reagents

- Sterile diluent for decimal dilutions: tryptone salt (TS)
- 16 \times 160 sterile tubes pre filled with 9ml of sterile TS
- Cooled crystal violet and neutral red lactose bile agar (VRBL).

10.4. Procedure

10.4.1. Agar medium

- Once prepared, keep the VRBL agar cooled in the water bath at $47^{\circ}\text{C} \pm 2^{\circ}\text{C}$ (for immediate usage).
- Never use a culture medium at a temperature higher than 50°C .
- Do not cool over 8 hours.
- For a deferred usage, keep agar cooled in an incubator at $55^{\circ}\text{C} \pm 1^{\circ}\text{C}$.
- Melted culture mediums unused within 8 hours, shall never re-solidify for later usage.

10.4.2. Culture

- Homogenise each dilution before inoculating in Petri dishes prior to carrying out decimal dilutions.
- Transfer 1 ml of reference solution and/or decimal dilutions retained in respective Petri dishes with pipettes changed after each dilution
- Introduce up to 20 minutes after the inoculum 15 to 20 ml of VRBL maintained in the water bath at $47^{\circ}\text{C} \pm 2^{\circ}\text{C}$.
- Slowly homogenise by shaking.
- Let solidify on laboratory bench (lid upwards).
- Pour approximately 5 ml of VRBL maintained in the water bath at $47^{\circ}\text{C} \pm 2^{\circ}\text{C}$.
- Let solidify on laboratory bench (lid upwards).
- Turn over dishes and incubate immediately 24 hours \square 2 hours at $30^{\circ}\text{C} \square 1^{\circ}\text{C}$.

10.4.3. Count

Dishes containing less than 150 characteristic or non-characteristic colonies based on two successive dilution are retained, but one of them must contain at least 15 characteristic colonies.

If only the dish inoculated with 1 ml of the 1st dilution contains under 15 characteristic colonies, then the dish will be retained for counting.

Characteristic colonies are counted manually or by using a counter.

Characteristic colonies after 24 hours \square 2 hours of incubation

- violet colonies surrounded sometimes by a red area (bile precipitation)
- diameter \approx 0.5 mm

10.5. Results

10.5.1. General case

Dishes containing less than 150 characteristic or non-characteristic colonies, based on two successive dilutions with one containing at least 15 characteristic colonies.

Method of calculation:

Number N of micro-organisms counted at 30 °C per millilitre (ml) or by gram (g) of product is obtained by calculating the weighted average of 2 retained dishes.

$$N = \frac{\sum c}{1,1d}$$

$\sum c$: sum of characteristic colonies counted of 2 retained dishes

d : dilution rate corresponding to the 1st dilution

Expression of results:

- round off the number N to the 2 largest digits
- -express to the tenth power

ex: $1.6 \cdot 10^3$ / g or ml

10.5.2. Estimation of small numbers

If the dish inoculated with 1 ml of the 1st dilution retained for analysis contains less than 15 characteristic colonies, the result will be expressed as follows:

$$N = c \frac{1}{d}$$

c : sum of characteristic colonies counted

d : rate of dilution

If the dish inoculated with 1ml of the 1st dilution retained for analysis contains no colonies then the result will be expressed as follows:

$$N = < 1 \frac{1}{d} \text{ micro-organisms per g or ml}$$

d : rate of dilution.

^[1] Point to be studied at a later date by the sub-commission of methods of analysis and appraisal of wine.

^[2] Point to be studied at a later date by the experts group “Wine microbiology”.