

**COMPENDIUM OF INTERNATIONAL ANALYSIS OF
METHODS – OIV
 α -dicarbonyl compounds by GC after derivatization**

**METHOD FOR THE DETERMINATION OF α -
DICARBONYL COMPOUNDS OF WINE BY GC AFTER
DERIVATIZATION BY 1,2-DIAMINOBENZENE
(OIV-Oeno 386B-2010)**

Method OIV-MA-AS315-21

Type IV method

1. Introduction

The principal α -dicarbonyl compounds found in wine (Fig 1) are: glyoxal, methylglyoxal, diacetyl and 2,3-pentanedione, but only α -diketones are relatively abundant in wine. Carbonyl compounds exist in all types of wines, particularly after malolactic fermentation and in red wines. In addition, sweet white wines produced with botrytized grapes can contain high levels of glyoxal and methylglyoxal.

Glyoxal: $\text{OCH}-\text{CHO}$ (ethanedial)

Methylglyoxal: $\text{CH}_3-\text{CO}-\text{CHO}$ (2-oxopropanal)

Diacetyl: $\text{CH}_3-\text{CO}-\text{CO}-\text{CH}_3$ (2,3-butanedione)

2,3-pentanedione: $\text{CH}_3-\text{CH}_2-\text{CO}-\text{CO}-\text{CH}_3$

2,3-hexanedione: $\text{CH}_3-\text{CH}_2-\text{CH}_2-\text{CO}-\text{CO}-\text{CH}_3$

Figure 1. The principal α -dicarbonyl compounds of wine (2,3-hexanedione is not naturally present in wine but it is used as internal standard).

**COMPENDIUM OF INTERNATIONAL ANALYSIS OF
METHODS – OIV
 α -dicarbonyl compounds by GC after derivatization**

Dicarbonyl compounds are important in wine for different reasons: their sensory impact, the reactivity with other components of the wine or possible microbiological effects.

2. Applicability

This method applies to all types of wines (white, red, sweetened or fortified), for carbonyl derivatives content ranging from 0.05 mg/L and 20 mg/L.

3. Principle

The method is based on the formation of derivatives of the quinoxaline type based on the α -dicarbonyl compounds of the wine with 1,2-diaminobenzene (Figure 2).

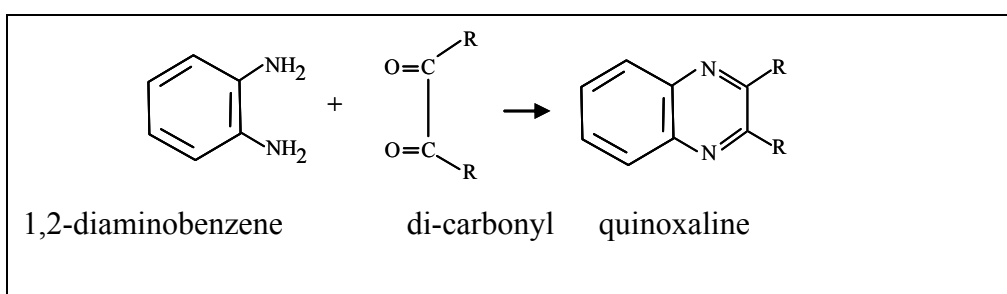


Figure 2 Formation of derivatives.

The reaction takes place directly in the wine at pH 8 and after a reaction time of 3 h at 60°C. The analysis of the derivatives is then carried out after extraction of the derivatives by dichloromethane and analysis by gas chromatography with detection by mass spectrometry (GC-MS) or using a nitrogen-specific detector.

4. Reagents and products

4.1 Dicarbonyl compounds

4.1.1 Glyoxal in a solution at 40% (CAS n° 107-22-3)

4.1.2 Methylglyoxal in a solution at 40% (CAS n° 78-98-8)

4.1.3 Diacetyl, purity > 99% (CAS n° 431-03-8)

4.1.4 2,3-Pentanedione, purity > 97% (CAS n° 600-14-6)

4.1.5 2,3-Hexanedione, purity > 90% (CAS n° 3848-24-6)

4.2 1,2-Diaminobenzene in powder form, purity > 97% (CAS n° 95-54-5)

**COMPENDIUM OF INTERNATIONAL ANALYSIS OF
METHODS – OIV
 α -dicarbonyl compounds by GC after derivatization**

4.3 Water for HPLC (for example microfiltered and with a resistivity of 18.2

M Ω)

4.4 Pure ethanol for HPLC (CAS n° 64-17-5)

4.5 Sodium hydroxide M. (CAS n° 1310-73-2)

4.6 Sulphuric acid 2M (CAS n° 7664-93-9)

4.7 Dichloromethane (CAS n° 75-09-2)

4.8 Anhydrous sodium sulphate (CAS n° 7757-82-6)

4.9 Aqueous-alcoholic solution at 50% vol .

Mix 50 ml of pure ethanol for HPLC (4.4) with 50 ml of water (4.3)

4.10 Solution of internal standard 2,3-hexanedione at 2.0 g/L

Place 40 mg of 2,3-hexanedione (4.2) in a 30-ml flask, dilute in 20 ml of aqueous-alcoholic solution to 50% vol (4.9) and stir until it has completely dissolved.

4.11 Anhydrous sodium sulphate (CAS n° 7757-82-6)

5. Equipment

5.1 Gas chromatograph with detection by mass spectrometry (GC-MS) or a nitrogen-specific detector.

5.1.1 Relatively polar, polyethylene glycol capillary column (CW 20M, BP21 etc.) with the following characteristics (as an example): 50 m x 0.32 mm x 0.25 μ m.

5.1.2 Data acquisition system.

5.2 pH measuring apparatus

5.3 Magnetic stirrer

5.4 Balance with a precision of 0.1 mg.

5.5 Oven which can be set to 60°C

5.6 Standard laboratory glassware including pipettes, screw-cap flasks, and microsyringes.

6. Preparation of the sample

No specific preparation is necessary

7. Procedure

Place 50 ml of wine in a flask (5.6)

**COMPENDIUM OF INTERNATIONAL ANALYSIS OF
METHODS – OIV
 α -dicarbonyl compounds by GC after derivatization**

Bring to pH 8 while stirring, with sodium hydroxide M (4.5)

Add 25 mg of 1,2-diaminobenzene (4.2)

Add 50 μ l of 2,3-hexanedione (internal standard) at 2.0 g/L (4.10)

Close the flask using a screw-cap fitted with a Teflon-faced seal

Stir until the reagent has completely disappeared (5.3)

Place in the oven at 60°C for 3 h (5.5)

Cool.

7.1 Optimisation and analytical conditions (this study was carried out by HPLC analysis, see this method)

The yield of the formation of derivatives of the dicarbonyl compounds with the 1-2-diaminobenzene is optimal at pH 8 at 60°C after three hours of reaction time

In addition, no interference of SO₂ with the formation of quinoxalines was noted during the study of the method.

7.2 Analysis by GC

7.2.1 Extraction of quinoxalines

- The reaction medium prepared in 7 is brought to pH 2 using H₂SO₄ 2M (4.6);
- Extract 2 times using 5 ml of dichloromethane (4.7) by magnetic stirring for 5 minutes;
- Decant the lower phase each time;
- Mix the two solvent phases;
- Dry on approximately 1 g of anhydrous sodium sulphate (4.11);
- Decant.

7.2.2 Chromatographic analysis (given as an example)

- *Detection.* For the analysis by GC-MS, a Hewlett Packard HP 5890 gas-phase chromatograph was coupled with Chemstation software and an HP 5970 mass spectrometer (electronic impact 70eV, 2.7 kV),

Note: It is also possible to use a nitrogen-specific detector

- *Column.* The column is a BP21 (SGE, 50 m x 0.32 mm x 0.25 μ m).

- *Temperatures.* The temperature of the injector and the detector are respectively 250°C and 280°C; that of the oven is held at 60°C for 1min,

**COMPENDIUM OF INTERNATIONAL ANALYSIS OF
METHODS – OIV
 α -dicarbonyl compounds by GC after derivatization**

then programmed to increase at a rate of 2°C/min to 220°C and the final isothermal period lasts 20 min.

- *Injection.* The volume injected is 2 μ l and the splitless time of the injector valves is 30s.

7.2.3 Analysis of quinoxalines formed

- *Separation.* The chromatogram of the derivatives from a wine obtained with 1,2-diaminobenzene, using selected-ion monitoring (SIM), is shown in Figure 3. Good separations were obtained with all types of wines (white, red, sweetened or fortified), and even with fermenting musts.

- *Identification of the peaks.* GC-MS was used to identify the dicarbonyl compounds derivatized from the wine based on the total ion current method (scan) which is used to obtain the mass spectra of derivatized quinoxalines and to compare them with those recorded in the library; in addition, the retention times were compared with those for pure compounds treated in the same way. Table 1 shows the principal ions of the mass spectra for the derivatized dicarbonyl compounds obtained.

- *Determination.* The quantitative determination of the dicarbonyl compounds is performed with the SIM method, by selecting ions $m/z = 76, 77, 103, 117, 130, 144, 158$ and 171 . The ions $m/z = 76$ and 77 are used for the quantification and the others as qualifiers, i.e. glyoxal: ions $m/z = 103$ and 130 , methylglyoxal: ions $m/z = 117$ and 144 , diacetyl: ions $m/z = 117$ and 158 , 2,3-pentandione: ions $m/z = 171$ and 2,3-hexanedione: ions $m/z = 158$ and 171 .

7.2.4 Characteristics of the method

Some elements of internal validation were determined, but this is not a formal validation according to the protocol governing the planning, the implementing and the interpreting of the performance studies pertaining to the analysis methods (OIV 6/2000)

- *Repeatability.* The repeatability of the GC-MS-SIM method shows coefficients of variation ranging between 2 and 5% for the four dicarbonyl compounds;

- *Recovery rate.* The quantities added to a wine were recovered with a recovery rate ranging between 92 and 117%;

**COMPENDIUM OF INTERNATIONAL ANALYSIS OF
METHODS – OIV
 α -dicarbonyl compounds by GC after derivatization**

- *Linearity*. Linear correlations were obtained in concentrations ranging from 0.05 to 20 mg/l.
- *Limit of detection*. The limit of detection of most of the derivatized dicarbonyl compounds using wine as a matrix is 0.05 mg/l

Table 1. Mass spectra (ion m/z and abundance of the ion in relation to that of the base peak) of derivatives of dicarbonyl compounds using 1,2-diaminobenzene

Dicarbonyl compound	Derivative	Mass spectrum (principal ions and abundance)
Glyoxal	Quinoxaline	130 (100), 103 (56.2), 76 (46.8), 50 (20.2), 75 (10.4), 131 (9.4)
Methylglyoxal	2-Methylquinoxaline	144 (100), 117 (77.8), 76 (40.5), 77 (23.3), 50 (21.9), 75 (11.3), 145 (10.3)
Diacetyl	2,3-Dimethylquinoxaline	117 (100), 158 (75.6), 76 (32.3), 77 (23.1), 50 (18.3), 75 (10.4)
2,3-Pentanedione	2-Ethyl-3-methylquinoxaline	171 (100), 172 (98), 130 (34.1), 75 (33.3), 77 (21), 50 (19.4), 144 (19), 143 (14.1), 103 (14)
2,3-Hexanedione	2,3-Diethylquinoxaline	158 (100), 171 (20.1), 76 (13.7), 77 (12.8),

**COMPENDIUM OF INTERNATIONAL ANALYSIS OF
METHODS – OIV**
 α -dicarbonyl compounds by GC after derivatization

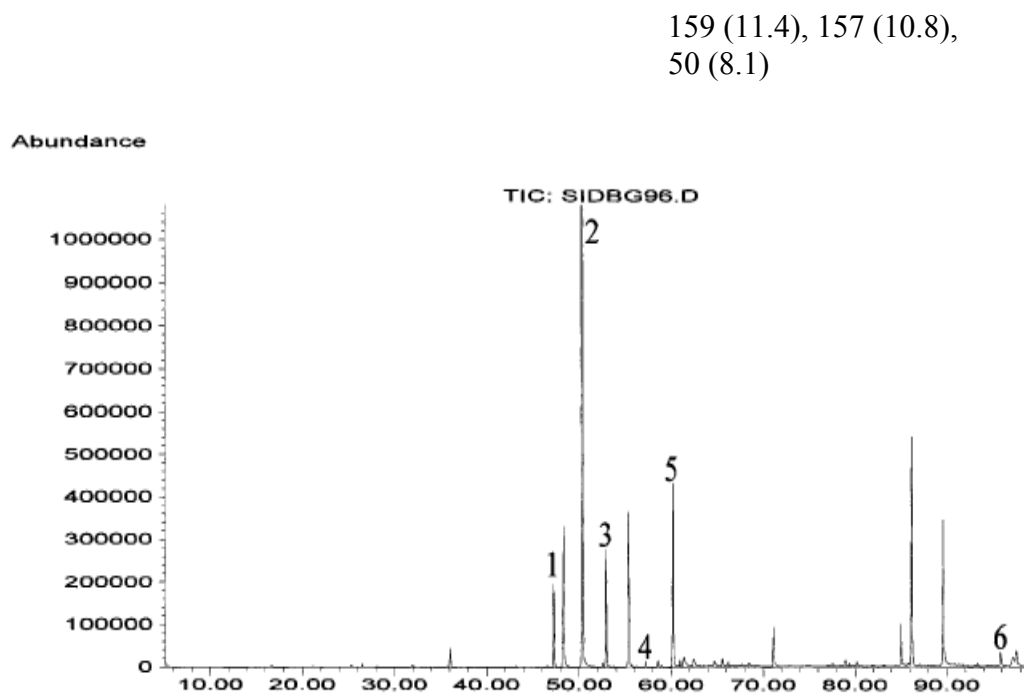


Figure 3. Gas chromatogram of the extract from the dicarbonyl compounds derivatized by 1,2-diaminobenzene from a white wine, detected by mass spectrometry by selecting the ions $m/z = 76, 77, 103, 117, 130, 131, 144, 158, 160$ and 171 . BP21 Column, $50\text{m} \times 0.32\text{mm} \times 0.25\text{ }\mu\text{m}$ oven temperature 60°C for 1min, then programmed increase of $2^\circ\text{C}/\text{min}$ up to 220°C . Injector temperature: 250°C .

1. glyoxal; 2. methylglyoxal; 3. diacetyl; 4. 2,3-pentanedione; 5. 2,3-hexanedione (internal standard); 6. phenylglyoxal (not studied with this method).

**COMPENDIUM OF INTERNATIONAL ANALYSIS OF
METHODS – OIV
 α -dicarbonyl compounds by GC after derivatization**

Bibliography

Bartowski E.J. and Henschke P.A. The buttery attribute of wine – diacetyl – desirability spoilage and beyond. *Int. J. Food Microbiol.* 96: 235-252 (2004).

Bednarski W., Jedrychowski L., Hammond E., and Nikolov L., A method for determination of α -dicarbonyl compounds. *J. Dairy Sci.* 72:2474-2477 (1989).

Leppanen O., Ronkainen P., Koivisto T. and Denslow J. A semiautomatic method for the gas chromatographic determination of vicinal diketones in alcoholic beverages. *J. Inst. Brew.* 85:278- 281 (1979).

Martineau B., Acree T. and Henick-Kling T., Effect of wine type on the detection threshold for diacetyl. *Food Res. Int.* 28:139-143 (1995).

Moree-Testa P. and Saint-Jalm Y., Determination of α -dicarbonyl compounds in cigarette smoke. *J. Chromatogr.* 217:197-208 (1981).

De Revel G., Pripis-Nicolau L., Barbe J.-C. and Bertrand A., The detection of α -dicarbonyl compounds in wine by the formation of quinoxaline derivatives. *J. Sci. Food Agric.* 80:102-108 (2000).

De Revel G. and Bertrand A. Dicarbonyl compounds and their reduction products in wine. Identification of wine aldehydes. *Proc 7th Weurman Flavour Research Symp., Zeist, June, pp 353-361 (1994).*

De Revel G. and Bertrand A., A method for the detection of carbonyl compounds in wine: glyoxal and methylglyoxal. *J. Sci. Food Agric.* 61:267-272 (1993).

Voulgaropoulos A., Soilis T. and Andricopoulos N., Fluorimetric determination of diacetyl in wines after condensation with 3,4-diaminoanisole. *Am. J. Enol. Vitic.* 42:73-75 (1991).

Gilles de Revel et Alain Bertrand, Analyse des composés α -dicarbonyles du vin après dérivation par le 1-2-diaminobenzène OIV FV 1275